

Simplify Your Data Review

waters_connect™ for Research, Investigations and Laboratories (RIL) and the QUAN Review Application: Advanced Tools for Managing LC-MS Data Analysis

Simplify Your LC-MS Data Processing

The QUAN Review application, part of the waters_connect for RIL product solution, provides a layered approach to data analysis in an easy-to-use package.



Key Features of QUAN Review:

- Reduce data review by up to 50%*
- The dashboard summarizes key parameters on one screen:
 - Standards, calibration, and QC are grouped together to streamline review.
- Exception Focused Review (XfR) allows users to easily see – and quickly address – errors that need attention.
- PanoGramic Display, a view that enables users to quickly identify problems or ambiguous data by showing multiple chromatograms, calibrations, and QC at the same time.

*As compared to TargetLynx™ Software

Set Flags or Alerts According to Your SOPs:

- Leverage the default result set to highlight data that is out of range such as R^2 , standard deviation, or signal-to-noise values.
- Select which rules are to be applied and the range/tolerance for each (e.g., the internal standard must be within 20% of the median response).
- Reduce turnaround time for labs with high sample volume and complex data sets.

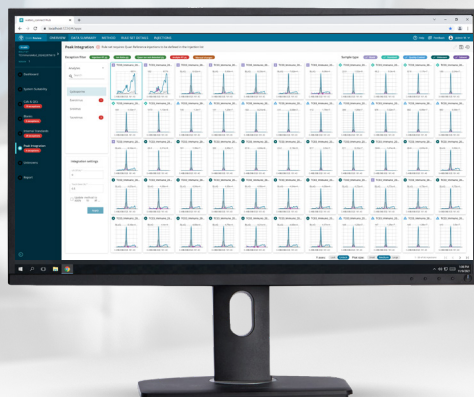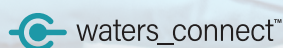

SCAN TO LEARN MORE

or contact us today at 360-567-8064 and unlock the full potential of your LC-MS/MS laboratory.

Waters™

Waters, ACQUITY, UPLC, waters_connect, and TargetLynx are trademarks of Waters Corporation. All other trademarks are the property of their respective owners.

©2024 Waters Corporation. June 24-12933 720008409EN Rev. A