


InertSearch™ for GC

InertCap® Applications

Basic samples

Data No. GA054-0000

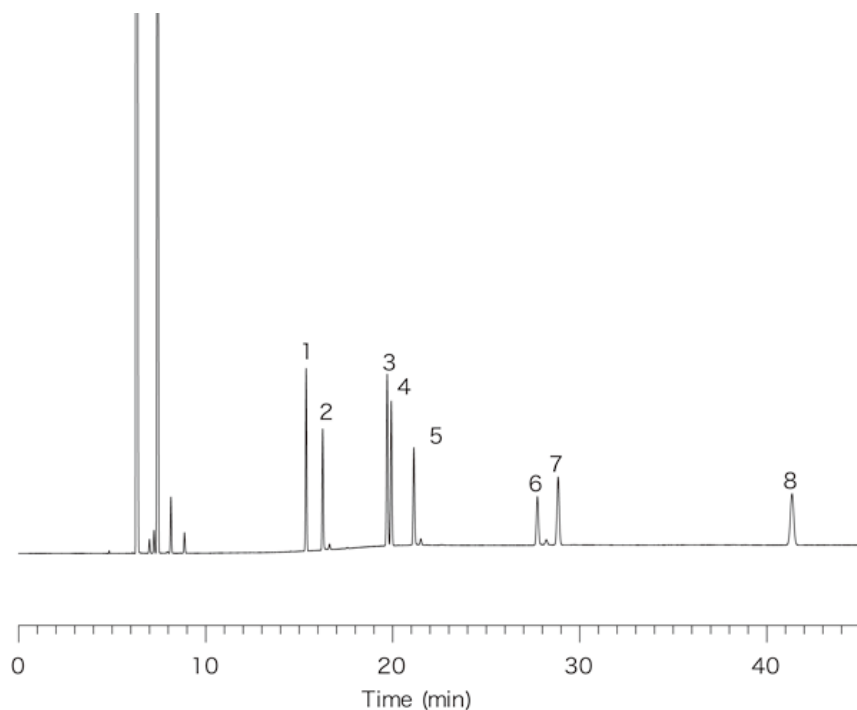

System :GC/FID
Column :InertCap for Amines
0.32mmI.D.×60m
Column Temp. :100°C→10°C/min→250°C (30min)
Carrier Gas :He 150kPa
Injection :Split 1 : 50
Detection :FID
Sample :Basic amines,Hydrocarbons

Solutes : 1. n-Decane
2. Amino octane
3. 2,6-Dimethylaniline
4. n-Dodecane
5. 1-Amino decane
6. Dicyclohexylamine
7. 1-Amino dodecane
8. n-Heptadecane