

GC CAPILLARY COLUMN APPLICATIONS: TRIGLYCERIDES

Italian Olive Oil

COLUMN: 007-65HT (65% Phenyl)
Methylpolysiloxane
15M. x 0.25mm I.D. x 0.1µm film

Cat. No.: 400-65HT-15-0.1F


Temperature: 350° (1.0 min. hold)
(0.5°/min.) - 365°C

Injector: 400°C

Detector: 400°C, FID

Carrier Gas: 35 cm/sec., Helium

1. pop
2. plp
3. pos
4. poo
5. plo
6. soo
7. ooo
8. olo
9. oll
10. aoo

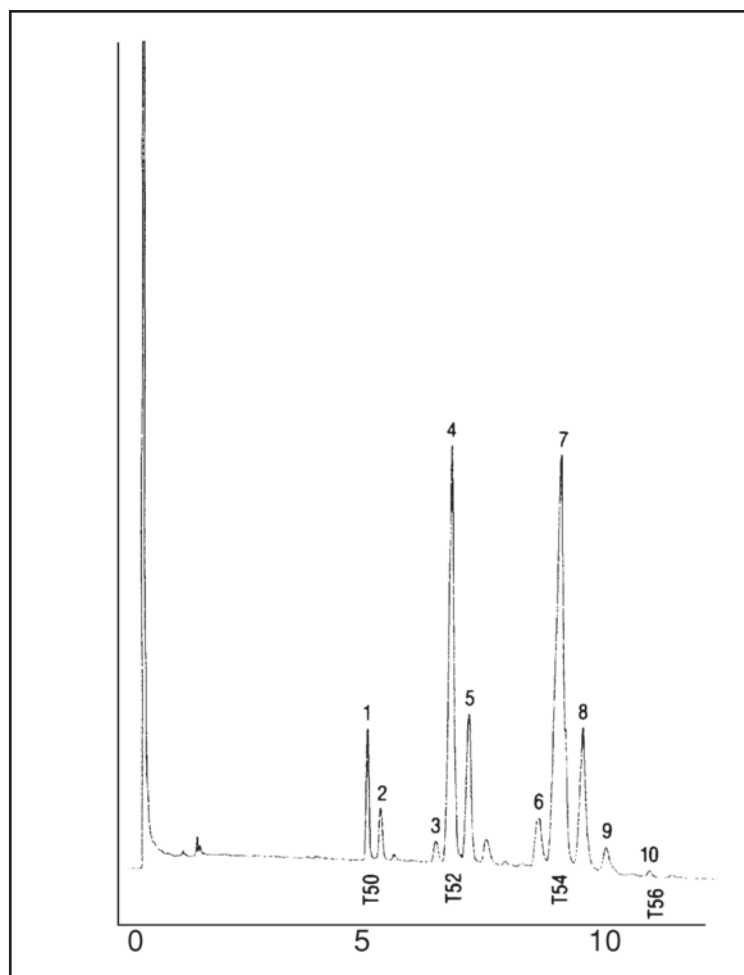

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com