



Analysis of Fluid Cat Cracker Feed using an Agilent J&W FactorFour VF-5ht UltiMetal Column

Application Note

Author

John Oostdijk
Agilent Technologies, Inc.

Introduction

This analysis of fluid cat cracker feed is performed using a VF-5ht UltiMetal column. The column has been developed using proprietary UltiMetal technology that provides a virtually unbreakable metal column material with excellent inertness properties similar to fused silica tubing. The UltiMetal tubing is coated with the VF-5 low bleed arylene stabilized liquid phase, resulting in a highly temperature stable and durable column perfectly suited for a variety of high temperature applications.



Agilent Technologies

Conditions

Technique:	GC-FID	Carrier Gas:	Hydrogen, constant flow mode	Injection Volume:	1.0 μ l
Column:	VF-5ht UltiMetal, 15 m x 0.32 mm Df = 0.1 μ m + Retention Gap, 2 m x 0.53 mm (p/n CP9095)	Injector:	On-column (1093), reversed liner, 100 $^{\circ}$ C (0 min) to 400 $^{\circ}$ C with 15 $^{\circ}$ C/min	Temperature:	50 $^{\circ}$ C (1 min) to 450 $^{\circ}$ C (20 mins) with 10 $^{\circ}$ C / min
Sample:	1 % CDU5 FCC Feed, 1 % in CS2	Detection:	FID (HT), 400 $^{\circ}$ C		

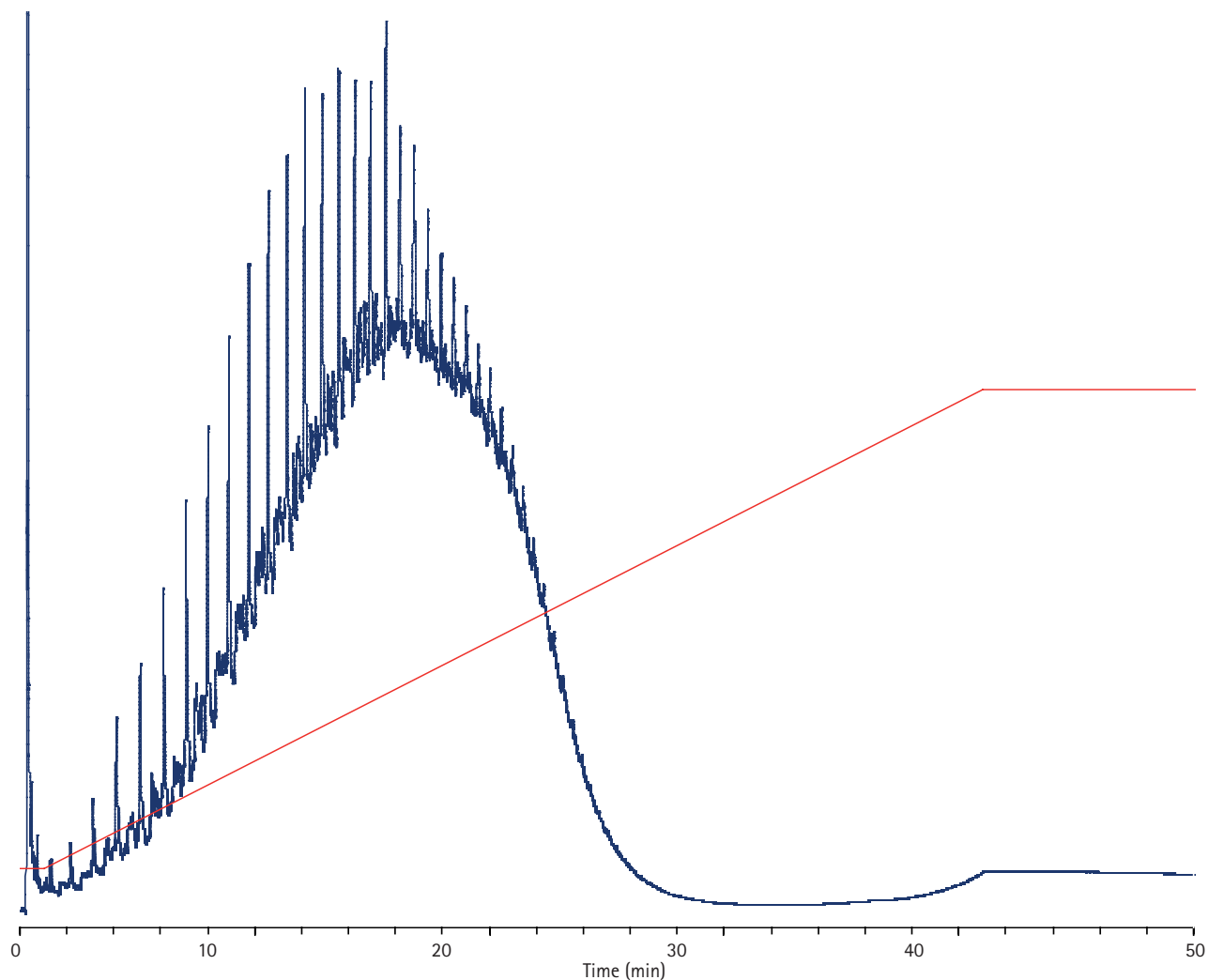


Figure 1. Analysis of fluid cat cracker feed using a VF-5ht UltiMetal column

www.agilent.com/chem

This information is subject to change without notice.

© Agilent Technologies, Inc. 2010

Published in UK, October 07, 2010

SI-01130



Agilent Technologies