

Separation of Alkyl Acetates

Method VEV0030J
HPLC RP Mode

Column: Eurospher II 100-3 C18, 150 x 3.0 mm ID

Order No. 10CE181E2G

Phase: Eurospher II 100-3 C18

Conditions: Eluent: Methanol / Water 60:40 (v/v)
Gradient: isocratic
Flow rate: 0.35 ml/min
Temperature: 25 °C
Volume: 5 µl


Detection: UV at 210 nm

Substances: Ethyl acetate, Butyl acetate, Pentyl acetate


Keywords: Alkyl acetates,

Chromatogram:


- 1 Ethyl acetate
- 2 Impurity A
- 3 Butyl acetate
- 4 Impurity B
- 5 Pentyl acetate

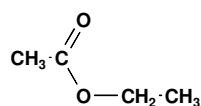

1 Ethyl acetate

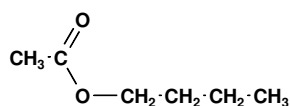

2 Butyl acetate

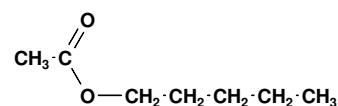

3 Pentyl acetate

